

METRICS

	∅ mm	G	↔ mm	PORT n.	n / f	REF			
INFILTRATION									
Skatter	1.5	17	100	5	∩	DOIN1510-2	10		
	2.0	14	150	5	∩	DOIN2015-2	10		
	2.5	13	250	5	∩	DOIN2525-2	10		
LIPOSUCTION									
Dolfin	3.0	11	150	3	∩	DOLF3150-2	10		
	3.0	11	250	3	∩	DOLF3250-2	10		
	4.0	8	150	3	∩	DOLF4150-2	10		
	4.0	8	250	3	∩	DOLF4250-2	10		
	4.0	8	350	3	∩	DOLF4350-2	10		
	5.0	6	250	3	∩	DOLF5250-2	10		
	5.0	6	350	3	∩	DOLF5350-2	10		
Dolfin Plus	4.0	8	250	7	∩	DOPL4250-2	10		
Dolfin Pro	4.0	8	250	8	∩	DOLP4250-2	10		
	4.0	8	350	8	∩	DOLP4350-2	10		
	5.0	6	250	8	∩	DOLP5250-2	10		
	5.0	6	350	8	∩	DOLP5350-2	10		
FAT GRAFTING									
HARVESTING									
Picking	MEGA	3.5	10	180	5	∩	PICK3518-2	10	
	Picking II	MEGA	4.0	8	300	8	∩	PICK352A-2	10
	MACRO	3.5	10	230	18	∩	PICK352L-2	10	
INJECTION									
Fillin Rigid	MEGA	3.5	10	250	1	∩	FIL1270S-2	10	
	MEGA	2.5	12	200	1	∩	FILU2515-2	10	
	MEGA	2.7	12	200	1	∩	FIL2590S-2	10	
	MEGA	2.7	12	150	1	∩	F9027150-2	10	
	MACRO	2.0	14	150	1	∩	FIL25C5S-2	10	
	MACRO	2.0	14	150	1	∩	FIL1670S-2	10	
	MACRO	2.0	14	100	1	∩	F9020100-2	10	
	MACRO	1.5	17	100	1	∩	F9015100-2	10	
	MACRO	1.5	17	70	1	∩	F9015070-2	10	
	MACRO	1.2	18	100	1	∩	F9012100-2	10	
	MICRO	1.2	18	70	1	∩	F9012070-2	10	
	Fillin Flexible	MICRO	1.0	19	100	1	∩	F9010100-2	10
		MICRO	1.0	19	70	1	∩	F9010070-2	10
MICRO		0.8	21	80	1	∩	F9008080-2	10	
MICRO		0.8	21	60	1	∩	F9008060-2	10	
MICRO		0.8	21	40	1	∩	F9008040-2	10	
NANO		0.6	23	50	1	∩	F9006050-2	10	
NANO		0.5	25	40	1	∩	F9005040-2	10	
NANO		0.5	25	25	1	∩	F9005025-2	10	
KITS									
		A	B1	B2	B3	C1	C2		
Macro St'rim	∅ ↔ mm	3.5 mm 220	12G 150	14G 150	17G 100	17G 70	PICS3523-2	5	
St'rim	∅ ↔ mm	2.0 mm 130	19G 50	21G 60	21G 40	14G 40 21G 40	PICS2130-2	5	
Mini St'rim	∅ ↔ mm	2.0 mm 90	19G 50	21G 60	21G 40	14G 40 21G 40	PICS2090-2	5	
Mini St'rim	∅ ↔ mm	2.0 mm 70	19G 50	21G 60	21G 40	14G 40 21G 40	PICS2070-2	5	
ACCESSORIES									
Doto	DETACHABLE HANDLE FOR MACRO ST' RIM						DOTO0155-2	10	

∩ EFFICIENCY


LOW
MEDIUM
HIGH
VERY HIGH

f EXTRUSION FORCE

LOW
VERY LOW

Be eSTRIM

PERFORM EXTREMELY SAFE & EFFICIENT FAT GRAFTING

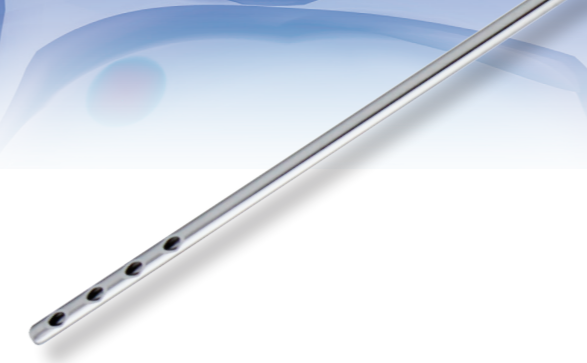

TYPE OF FAT GRAFTING	SIZING HARVESTING CANNULAS Picking's			EQUALIZING INJECTION CANNULAS Fillin			KITS HARVESTING & INJECTION
	Ref.	Lobule Size	∩	∅ mm	f	Ref.	
MEGAFAT	Picking	 < 2700 x 3800 μ	HIGH	10 G	3,5	VERY LOW	—
	Picking II		VERY HIGH	12 G	2,7	LOW	
MACROFAT	Picking II	 < 1000 x 1600 μ	HIGH	14G	2,0	VERY LOW	Macro St'rim
				17G	1,5	LOW	
MICROFAT		 < 450 x 600 μ		18G	1,2	VERY LOW	St'rim Mini St'rim
				19G	1,0	VERY LOW	
NANOFAT MAT ¹ /tSVF ²		 < 80 x 100 μ		21G	0,8	LOW	—
				23G	0,6	VERY LOW	
				25G	0,5	LOW	

LIPU_FG_PORTFOLIO.CPC.EN.C4(OCTOBER23)

The spheres are calibrated and not deformable and are only intended to illustrate a hypothesis of an adipocyte lobule size based on data of the technical files of STRIM HEALTH CARE.

- (1) Microfragmented Adipose Tissue
- (2) Stromal Vascular Fraction (tissue)

∩ Suction efficiency
f FILLIN'S extrusion forces

FIND OUT MORE ABOUT IN
WWW.STRIMHEALTHCARE.COM




welcome@strimhc.com | +33 450 72 17 18 | Boulevard des Philosophes 17 - Geneva -Switzerland

©2023

strim
HEALTH CARE

LIPOSUCTION & FAT GRAFTING


COMPACT PRODUCT CATALOG

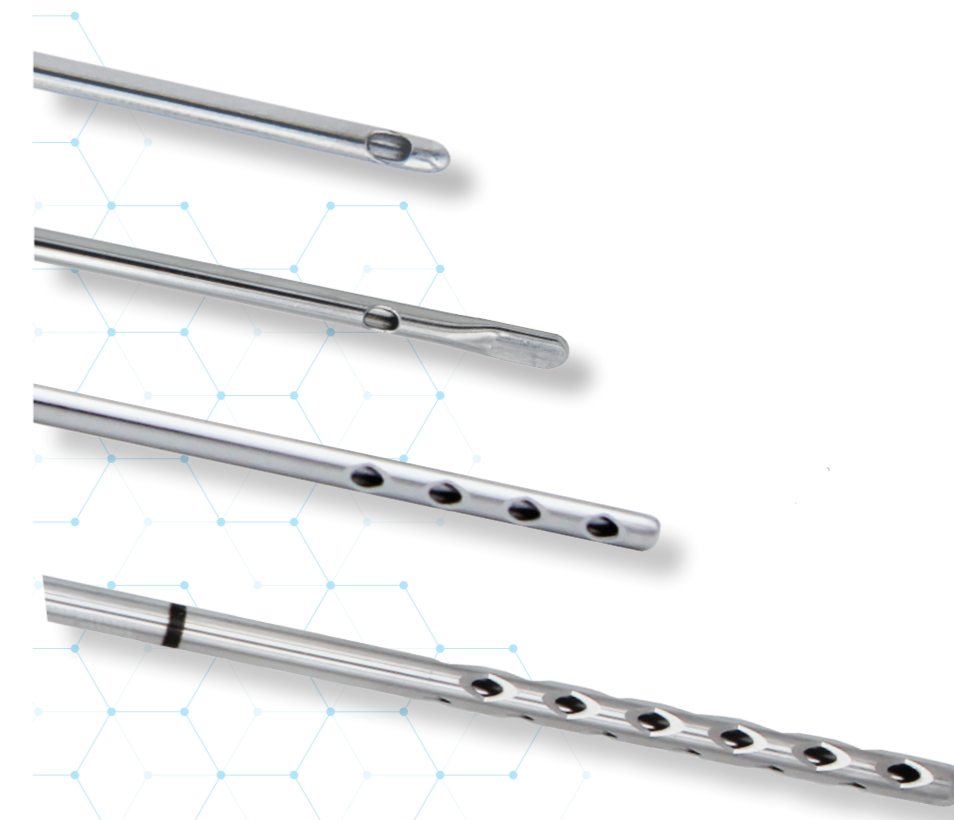


INFILTRATION

SKATTER'S


Multi-spiraled cannulas for tumescent fluid infiltration prior to fat harvesting

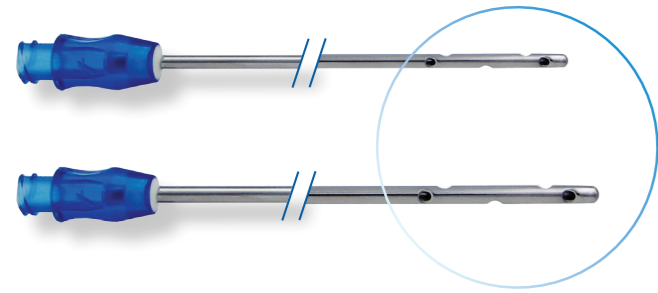

Length: 100, 150 and 250 mm

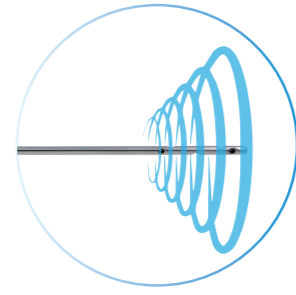

Vortex Effect

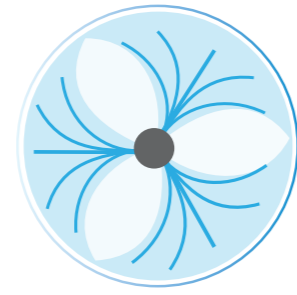

360° Distribution

LIPOSUCTION

DOLFIN'S

Body reshaping cannulas with sharp edges ports "guillotine" design to cut adipose tissue




DOLFIN


3 PORTS

1+1+1

MEDIUM η

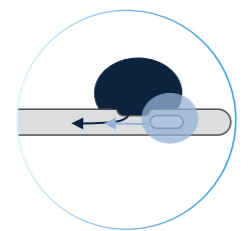

260° SUCTION PRECISION

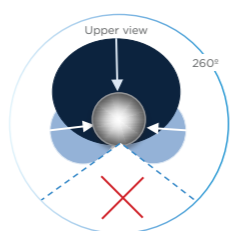


DOLFIN PLUS


7 PORTS

2+3+2

HIGH η

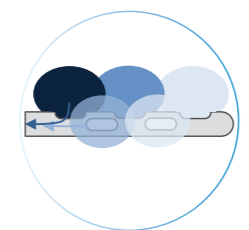

260° SUCTION PRECISION

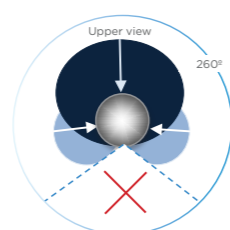


DOLFIN PRO


8 PORTS

4+4

VERY HIGH η

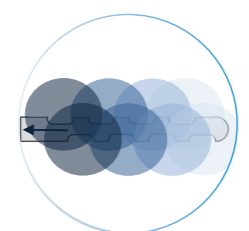

360° FULL SUCTION

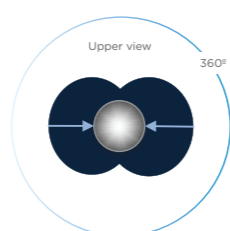

FAT GRAFTING

MEGAFAT



PICKING II
PICK352A-2

Ø mm	4.0	3.5	10G	12G	
← (mm)	300	180	250	200	150
PORTS	8	5	INJECTION		
η	VERY HIGH	HIGH			
SUCTION	260°	360°			

HARVESTING

HARVESTING

MACROFAT



PICKING II
PICK352L-2

Ø mm	3.5	14G	17G	
← (mm)	230	150	100	70
PORTS	18	INJECTION		
η	HIGH			
SUCTION	260°			

HARVESTING

ALL IN ONE

KITS

macro st'rim® MACROFAT GRAFTING KIT




OPTIONAL:

- Connector FLL/FLL
- DOTO handle
- SKATTER



MICROFAT

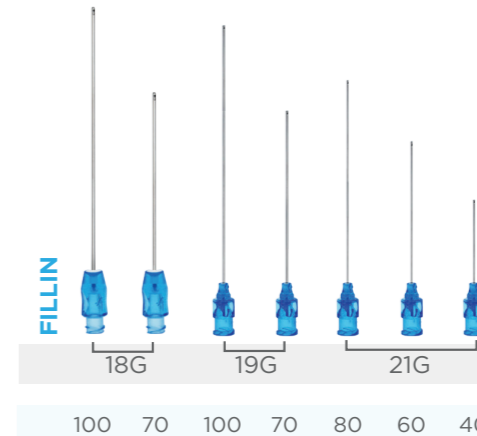


FILLIN

18G 19G 21G

100 70 100 70 80 60 40

NANOFAT

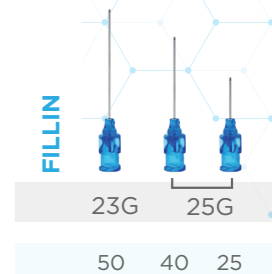


FILLIN

23G 25G

50 40 25

st'rim®

MICROFAT GRAFTING KIT

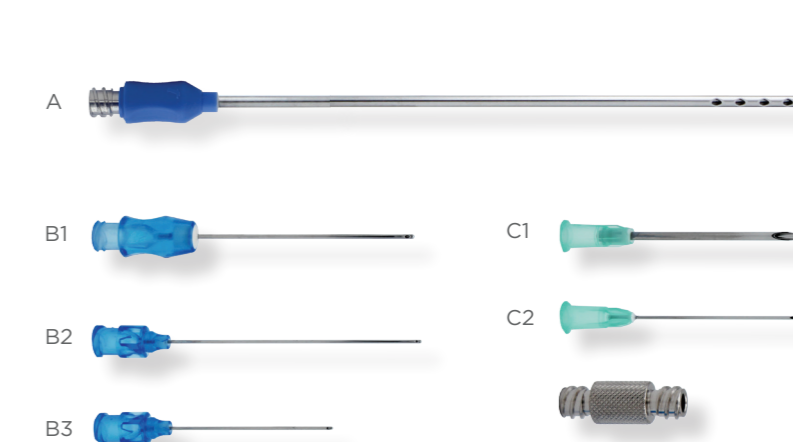

mini st'rim®

Similar kit to st'rim but with shorter harvesting cannulas - see metrics tab.